

Laufzeittabelle GAMMAT® 11C

Desilliertes Wasser (ml)
als Wärmebad im 170 ml Glas

Goldanteil (g)


ECOLYT SG100 (ml)

ACTIVATOR SG100 (ml)

Laufzeit (h)

70 ml Becherglas Rührer: Ø 7 x 25	60				
		0,60	40	0,50	7
		0,75	50	0,60	7
		0,90	60	0,75	7
		1,05	70	0,90	7
		1,20	80	1,00	7

170 ml Becherglas Rührer: Ø 6 x 35				
	1,50	100	1,40	8
	1,65	110	1,50	8
	1,80	120	1,70	8
	1,95	130	1,80	8
	2,10	140	2,00	8
	2,25	150	2,20	8
	2,40	160	2,30	8
	2,55	170	2,50	8
	2,70	180	2,70	8
	2,85	190	2,80	8
	3,00	200	3,00	8

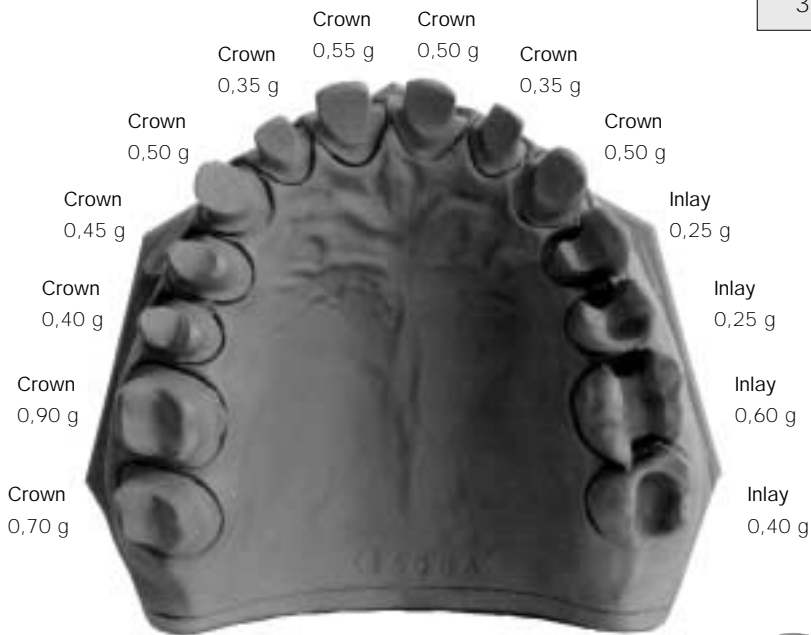

Brücke 3-gliedrig
1,80 g
100 ml



Brücke 3-gliedrig
2,10 g
120 ml

Gramm Technik
Dentale Galvanotechnik



Gramm GmbH
Postfach 1124
D-75229 Tiefenbronn-Mühlhausen
Telefon 07234/9519-0
Telefax 07234/9519-40
e-mail gramm.dental@t-online.de
<http://www.galvanoforming.de>